

Optima^{INT} Buying Guide



OPTIMA^{INT} MX2 PRO

The Optima^{INT} motor system simplifies restorative and endodontic procedures. Take control with a fully integrated, compact system that provides up to two MX2 PLUS LED micromotors, programmable presets, and WaveOne[®] Gold, compatible with RECIPROC[®]*, for simplified endodontics.

RESTORATIVE SIMPLIFIED

The MX2 PLUS LED micromotor delivers constant torque across its speed range of 100-40,000 RPM. Use a single 1:5 contra-angle attachment for nearly every restorative procedure for a full speed range of 500-200,000 RPM.

ROTARY ENDO WITH CONFIDENCE

Set the maximum allowed torque based on file manufacturer guidelines, and the system will auto-reverse the file upon exceeding the preset torque limits, preserving the file.

WAVEONE[®] GOLD & RECIPROC[®] READY


New to the DMX3 control board: WaveOne[®] Gold and RECIPROC[®]-ready functionality when paired with the CA ENDO 4:1 attachment (sold separately, REF 1600955-001).

REDUCE DOWNTIME: OPTIMA^{INT} MX2 PRO DUAL MOTOR KIT

Use two attachments of different speeds for specialized restorative procedures, or seamlessly switch between restorative and endo treatments in the same operation. Lift an attachment from the delivery unit to instantly recall its previous settings.

RESTORATIVE ONLY: OPTIMA^{INT} MCX SINGLE MOTOR KIT

For dentists who need restorative mode only, Bien-Air offers a unit with the MCX LED micromotor that provides equivalent performance with the same full integration at a lower price.

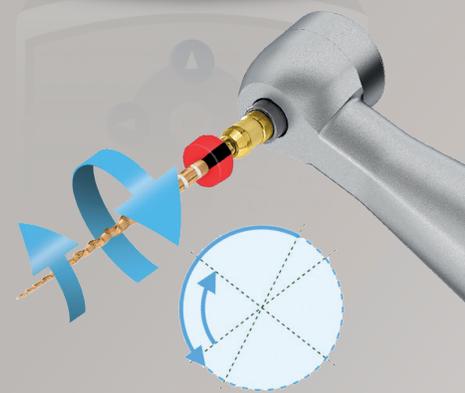

OPTIMA ^{INT} MCX SINGLE MOTOR KIT	REF 1700636-001
OPTIMA ^{INT} MX2 PRO SINGLE MOTOR KIT	REF 1700770-001
OPTIMA ^{INT} MX2 PRO DUAL MOTOR KIT	REF 1700771-001

*For RECIPROC[®] and RECIPROC[®] blue endo files, use the WaveOne[®] setting with CA ENDO (REF 1600955-001).

S W I S S  M A D E

Optima^{INT} Buying Guide



	1700636-001 Single MCX Motor	1700770-001 Single MX2 PLUS LED Motor	1700771-001 Dual MX2 PLUS LED Motor
Internal Installation	✓	✓	✓
Restorative	✓	✓	✓
Rotary Endodontics WaveOne® Gold/RECIPROC® Ready	✗	✓	✓

MOUNTING BRACKETS

Each Optima^{INT} motor system includes a bracket for touchpad installation onto most delivery units. For full integration into the following delivery units, order their corresponding brackets.

INCLUDED BRACKET (RIGHT SIDE MOUNT)	REF 1302927-001A
LEFT SIDE VERSION OF INCLUDED BRACKET	REF 1302927-001AL
ADEC 300	REF 1304437-001 & 1304437-HRD
ADEC 500	REF 1304437-001
ADEC 532B/533B (LEFT SIDE MOUNT ONLY)	REF 1308204-001
ADEC 541 REAR DELIVERY UNIT	REF 1900345-001

For delivery units made by Belmont, Midmark, and Ritter, contact the manufacturer for bracketing information. If you have a legacy delivery unit or are unsure, contact your local Bien-Air sales representative for additional custom bracket information.



Connect with us and schedule
a **FREE** in-office demo!
bienairusa.com/connect-with-us



WaveOne Gold is a registered trademark of DENTSPLY SIRONA. RECIPROC is a registered trademark of VDW.

S W I S S  M A D E